


SI-LUTION[®] SACROILIAC JOINT FUSION SYSTEM



Features & Benefits

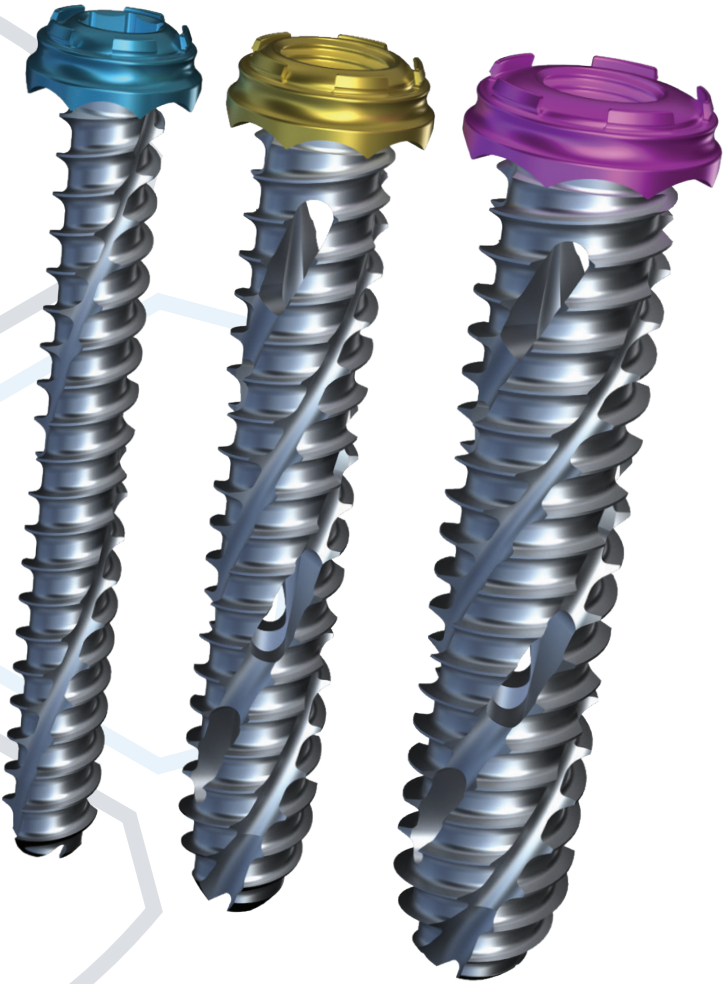

CURITEVA[®]



SI-LUTION[®] Sacroiliac Joint Fusion System

The SI-LUTION Sacroiliac Joint Fusion System is designed to offer immediate stability and long-term fixation across the sacroiliac (SI) joint. The implants are machined out of Titanium and incorporate the proprietary TRU-LOK[®] micro-textured surface technology which is intended to provide an improved osteogenic environment and encourage enhanced osseointegration. The system offers screws in multiple diameters and various lengths as well as stabilization caps to provide an ideal solution for every patient.

System Features

- Low profile head
- Partially fluted head to aid screw seating
- TRU-LOK[®] surface technology
- Optional color-coded stabilization caps provide increased load distribution

Stabilization Caps (Optional)

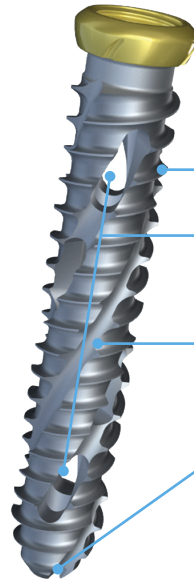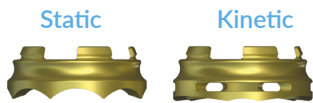

- Dual-lead threads for rapid screw advancement
- Large centralized helical fenestrations allow for direct allograft integration across the SI joint
- Self-harvesting bone channels to harvest and distribute local bone during insertion
- Aggressively tapered tip to facilitate implant insertion

SACROILIAC SCREWS		
Catalog Number	Description	Quantity
C614-7030	Fully Threaded Screw, 7 X 30mm	*
C614-7035	Fully Threaded Screw, 7 X 35mm	2
C614-7040	Fully Threaded Screw, 7 X 40mm	2
C614-7045	Fully Threaded Screw, 7 X 45mm	2
C614-7050	Fully Threaded Screw, 7 X 50mm	2
C614-7055	Fully Threaded Screw, 7 X 55mm	2
C614-7060	Fully Threaded Screw, 7 X 60mm	2
C614-7065	Fully Threaded Screw, 7 X 65mm	2
C614-7070	Fully Threaded Screw, 7 X 70mm	*
C614-9530	Fully Threaded Screw, 9.5 X 30mm	*
C614-9535	Fully Threaded Screw, 9.5 X 35mm	3
C614-9540	Fully Threaded Screw, 9.5 X 40mm	3
C614-9545	Fully Threaded Screw, 9.5 X 45mm	3
C614-9550	Fully Threaded Screw, 9.5 X 50mm	3
C614-9555	Fully Threaded Screw, 9.5 X 55mm	3
C614-9560	Fully Threaded Screw, 9.5 X 60mm	3
C614-9565	Fully Threaded Screw, 9.5 X 65mm	3
C614-9570	Fully Threaded Screw, 9.5 X 70mm	*

SACROILIAC SCREWS		
Catalog Number	Description	Quantity
C614-1230	Fully Threaded Screw, 7 X 30mm	*
C614-1235	Fully Threaded Screw, 7 X 35mm	2
C614-1240	Fully Threaded Screw, 7 X 40mm	2
C614-1245	Fully Threaded Screw, 7 X 45mm	2
C614-1250	Fully Threaded Screw, 7 X 50mm	2
C614-1255	Fully Threaded Screw, 7 X 55mm	2
C614-1260	Fully Threaded Screw, 7 X 60mm	2
C614-1265	Fully Threaded Screw, 7 X 65mm	2
C614-1270	Fully Threaded Screw, 7 X 70mm	*

STABILIZATION CAPS (OPTIONAL)		
Catalog Number	Description	Quantity
C614-70-SC	Static Stabilization Cap, 7mm	4
C614-95-SC	Static Stabilization Cap, 9.5mm	4
C614-12-SC	Static Stabilization Cap, 12mm	4
C614-70-KC	Kinetic Stabilization Cap, 7mm	*
C614-95-KC	Kinetic Stabilization Cap, 9.5mm	*
C614-12-KC	Kinetic Stabilization Cap, 12mm	*

*Available by request



For more information or to place an order call: 877.9CURITEVA, email: customersupport@curiteva.com or visit www.curiteva.com

Products patented and/or patent pending.
All products are not currently available in all markets.

Copyright©2022 by Curiteva, Inc.
QMSF-1187 Rev 2.0