


SAVANT™ PEEK

Posterior Lumbar Interbody Fusion System
Featuring Titanium Plasma Coating

Features & Benefits

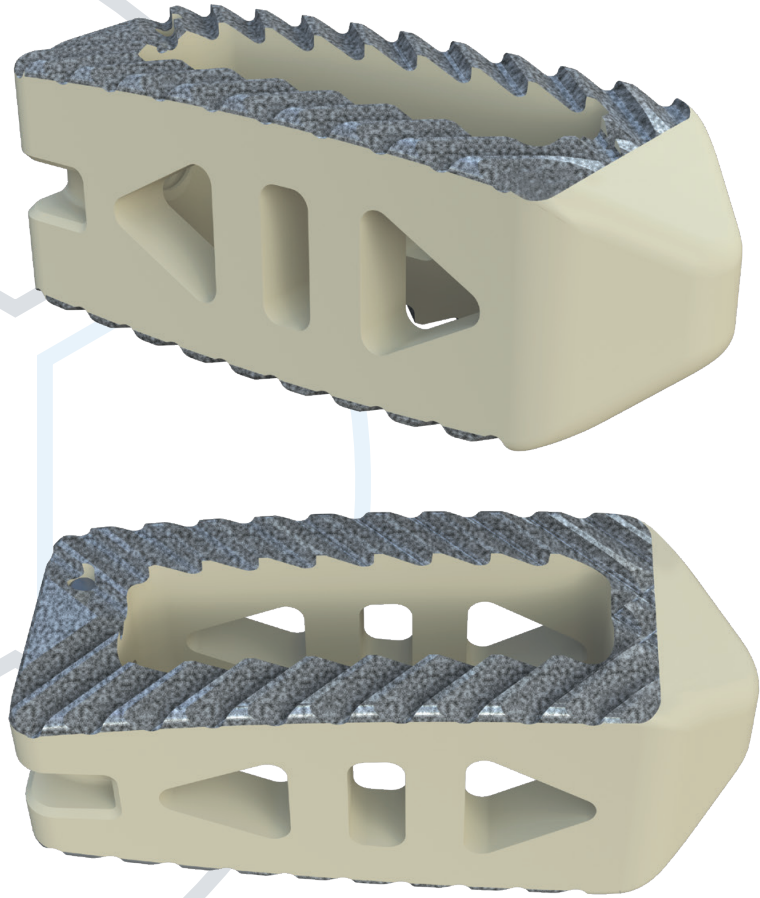


CURITEVA®

SAVANT™ PEEK PLIF System

Featuring Titanium Plasma Coating

The Savant Titanium Plasma Coated Posterior Lumbar Interbody Fusion System is designed to offer immediate stability and long-term fixation. The interbody is manufactured out of PEEK with a titanium plasma spray coating. The titanium porous coating creates topographical properties and provides fluoroscopic endplate visualization. The interbody is offered in multiple footprints, heights and lordosis providing an ideal fit for every patient.

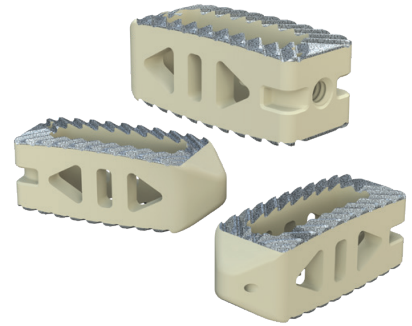

System Features:

- Surface topography produces superior implant contact and construct stability
- Modulus of elasticity between cortical and cancellous bone for optimal load sharing
- Contoured, distracting nose to facilitate placement
- Designed to ensure a high degree of compressive strength and dimensional stability
- Central and lateral graft windows for biologic seeding
- Multi-directional teeth resist migration
- Optimal Radiolucency
- Tantalum radiopaque markers to optimize visibility and placement
- Comprehensive discectomy instruments
- Full range of Trials for proper sizing of implant height

Titanium Plasma Coated PEEK PLIF Interbodies

22mm Length, 10mm Width, 0°		26mm Length, 10mm Width, 0°		30mm Length, 10mm Width, 0°	
Height	Catalog Number	Height	Catalog Number	Height	Catalog Number
6mm	C324-221006-0CT	6mm	C324-261006-0CT	6mm	C324-301006-0CT
7mm	C324-221007-0CT	7mm	C324-261007-0CT	7mm	C324-301007-0CT
8mm	C324-221008-0CT	8mm	C324-261008-0CT	8mm	C324-301008-0CT
9mm	C324-221009-0CT	9mm	C324-261009-0CT	9mm	C324-301009-0CT
10mm	C324-221010-0CT	10mm	C324-261010-0CT	10mm	C324-301010-0CT
11mm	C324-221011-0CT	11mm	C324-261011-0CT	11mm	C324-301011-0CT
12mm	C324-221012-0CT	12mm	C324-261012-0CT	12mm	C324-301012-0CT
13mm	C324-221013-0CT	13mm	C324-261013-0CT	13mm	C324-301013-0CT
14mm	C324-221014-0CT	14mm	C324-261014-0CT	14mm	C324-301014-0CT
15mm	C324-221015-0CT	15mm	C324-261015-0CT	15mm	C324-301015-0CT
16mm	C324-221016-0CT	16mm	C324-261016-0CT	16mm	C324-301016-0CT

22mm Length, 10mm Width, 8°		26mm Length, 10mm Width, 8°		30mm Length, 10mm Width, 8°	
Height	Catalog Number	Height	Catalog Number	Height	Catalog Number
7mm	C324-221007-8CT	8mm	C324-261008-8CT	9mm	C324-301009-8CT
8mm	C324-221008-8CT	9mm	C324-261009-8CT	10mm	C324-301010-8CT
9mm	C324-221009-8CT	10mm	C324-261010-8CT	11mm	C324-301011-8CT
10mm	C324-221010-8CT	11mm	C324-261011-8CT	12mm	C324-301012-8CT
11mm	C324-221011-8CT	12mm	C324-261012-8CT	13mm	C324-301013-8CT
12mm	C324-221012-8CT	13mm	C324-261013-8CT	14mm	C324-301014-8CT
13mm	C324-221013-8CT	14mm	C324-261014-8CT	15mm	C324-301015-8CT
14mm	C324-221014-8CT	15mm	C324-261015-8CT	16mm	C324-301016-8CT
15mm	C324-221015-8CT	16mm	C324-261016-8CT		
16mm	C324-221016-8CT				



For more information or to place an order call: 877.9CURITEVA, email: customersupport@curiteva.com or visit www.curiteva.com

Products patented and/or patent pending.
All products are not currently available in all markets.

Copyright©2020 by Curiteva, Inc.
QMSF-1147 Rev 1.0